

PROSAT® Sterile™ PS-911EB

- Nonwoven, meltblown polypropylene presaturated wipes
- Saturated with 70% isopropyl alcohol/30% deionized water
- Sterilized by Gamma Irradiation
- Validated sterile to a Sterility Assurance Level of 10^{-6} in accordance with AAMI Guidelines.
- Certificate of Sterility is available upon request
- Tamper evident label
- Resealable pouch preserves cleanliness and solvent saturation levels
- Each package is individually lot-numbered



Applications

- Designed for sterile environments
- ISO Class 5 compatible
- Routine cleaning and wipe down of lab tools, instruments and other equipment
- Cleanup following Rhodac plate sampling


Catalogs: [Products for Cleanrooms and Controlled Environments](#)

Data Sheets: [PS_911EB](#)

Part #	Description	Packaging	SDS	Sample
PS-911EB	9" x 11" (22.9 x 27.9cm) wipe, sterile	30 wipes/pouch, 12 pouches/box, 4 boxes/case; 48 pouches/case		

PROSAT® Sterile™ Polypropylene Low Endotoxin Wipes

Contec's PROSAT® Sterile Polypropylene Low Endotoxin wipes are presaturated with the optimal amount of 70% USP Grade IPA and 30% Deionized Water. The nonwoven meltblown polypropylene wipes provide a consistent release of solvent to thoroughly remove surface contaminants in critical environments. Contec's careful ingredient selection and strict process control ensure the wipes meet USP standards for low endotoxin products. Each lot is tested before release and low endotoxin certified to <20 EU/wipe, eliminating the risk of introducing endotoxins and other contaminants into product contact areas.

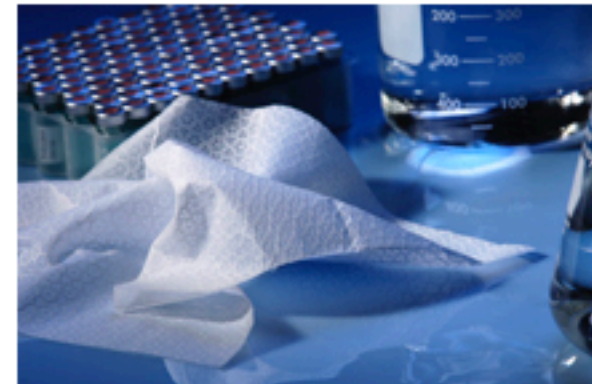


The resealable pouch reduces waste. Terminally sterilized by gamma irradiation, the wipes are Validated Sterile to a 10⁻⁶ Sterility Assurance Level.

Applications

- ISO Class 5 (Grades A/B) environments compatible
- Ideal for use in aseptic processing areas.

Part #	Description	Packaging	SDS	Sample
PS-911LE	 PROSAT Sterile Polypropylene LE Wipes, 9" x 11" (23 x 28cm), double bagged, sterile	30 wipes/pouch; 36 pouches/case		

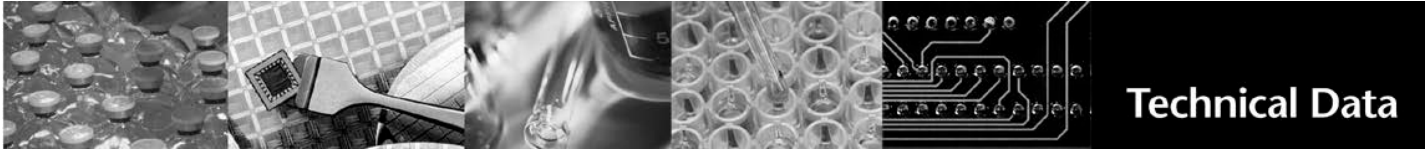

PROSAT® Sterile™ Wipes, Meltblown Polypropylene

Meltblown Polypropylene with 70% USP Grade IPA, 30% DI Water

Description

PROSAT *Sterile* Wipes consist of 100% polypropylene wipes saturated with a solution of 70% isopropanol alcohol and 30% deionized water in a resealable pouch. The pouch opens easily and the wipes are readily dispensed.

PROSAT *Sterile* Wipes utilize a meltblown polypropylene wipe which is exceptionally clean, being free from additives of any kind, and has a very low concentration of sodium. Increase solvent control and accountability and reduce VOC emissions with this sterile wipe. Compatible with ISO Class 5 cleanrooms.

Validated Sterile to a Sterility Assurance Level of 10^{-6} in accordance with AAMI guidelines.

Technical Data

Wipe material	Polypropylene	
Wipe construction	Meltblown nonwoven	
Solution of saturation; (blend ratio and constituents)	70% isopropanol alcohol/ 30% deionized water	
Solvent Grade	USP	
Attribute; (units)	Value	Test Method
Basis weight; (g/m ²)	36	Contec Method
Non-volatile residue, NVR		IEST-RP-CC004.3, Sec. 7.1.2
In deionized water; (g/m ²)	0.01	
In isopropanol; (g/m ²)	0.07	
Specific ions		IEST-RP-CC004.3, Sec. 7.2.2
Sodium; (ppm)	3.0	
Chloride; (ppm)	5.0	
Particles, readily releasable		IEST-RP-CC004.2, Sec. 5.1
P ≥ 0.5µm; (x10 ⁶ /m ²)	5.8	
Fibers > 100µm; (x10 ³ /m ²)	6.3	IEST-RP-CC004.2, Sec. 5.2
Saturation volume tolerance;		
< 300 mL	+/- 12%	
300 -1666 mL	+/- 3%	
> 1666 mL	+/- 50 mL	

Part Number	Description	EA/PCH	PCH/OB1	OB1/BX	BX/CS	EA/CS
PS-911EB	PROSAT <i>Sterile</i> wipes, 9" x 11" (23 x 28cm), 70% IPA/30% DI Water	30	6	2	4	1440
PSP0039	PROSAT <i>Sterile</i> wipes, 9" x 11" (23 x 28cm), 70% IPA/30% DI Water, updated configuration as of February 2016	50	4	1	6	1200
PS-840IR	PROSAT <i>Sterile</i> wipes, 8" x 8" (20 x 20cm), 70% IPA/30% DI Water	40	4	1	6	960

Key: EA = Each, OB1 = Outer Bag 1, PCH = Pouch, BX = Box, CS = Case

Notes:

- 1) The information presented here is applicable to the part numbers shown above as well as to any product containing the same materials and produced under the same conditions, regardless of product size or packaging configuration. Please contact a Contec sales representative for more details.
- 2) Data shown are typical values and should not be used as product specifications.
- 3) Valid product comparisons may only be obtained through side-by-side testing in the same test facility, under similar conditions.
- 4) Current and/or comparison data may be available. Please contact a Contec sales representative for details.

Revision date: 02/08/16



Copyright © 2015 Contec, Inc. All rights reserved.

North America
P.O. Box 530
Spartanburg, SC 29304 USA
Phone: +1-864-503-8333
www.contecinc.com

Europe
Zl du Prat – RP 3707
56037 VANNES Cedex
FRANCE
+33 297 437 690

China
No.17 Longyun Road
Suzhou 215024
CHINA
+ 86-512-6274 4050